

## GYNECOLOGISTS' BENEFITS

All gynecologists can perform dilation using Gynedil® with less risk and its applicability is independent of the provider's experience.

Gynedil is intended for single use, so there is reduced risk of intrauterine infection.

Only one insertion is performed and no redilatation is required.

Total comfort for the doctor because the risk of cervical tissue damage is significantly reduced.

Dilation with Gynedil takes between 20 to 30 seconds which shortens the duration of the whole procedure.

## PATIENTS' BENEFITS

Symptoms associated with "cervical incompetence" were found in 75% of women who had violent dilatation during termination of pregnancies<sup>1</sup>.

- Gynedil's build reflects an innovative approach that enables a more effective and less traumatic dilation.
- Its flexible balloon dilator allows penetration in the cervical canal with minimal resistance and adaptation to curvature of the patient's anatomy. This prevents the risk of infection, laceration of the cervix, perforation of the uterus and cervical incompetence.

1. A survey by Dr. Ren, "Cervical Incompetence - Aetiology and Management," Medical Journal of Australia (December 1993,29), Volume 60

+49 211 650 415 60

info@gynedil.com

www.gynedil.com

D.Med Healthcare  
Germany AG  
Klaus-Bungert-Straße 3  
40468 Düsseldorf  
Germany

# GYNEDIL®

## Hydraulic Cervical Dilator



### For Gynecologists:

An optimized medical device that gynecologists can rely on.


### For Women:

An innovative solution to explore towards making informed decisions on the best approach for cervical dilation.

A PRODUCT OF  
**D.med**  
healthcare

## Gynedil - Hydraulic dilator for the cervical canal

Gynedil consists of a manual hydraulic pump filled with distilled water and a flexible balloon dilator. The balloon dilator is 6cm long with an initial width diameter of 4.5mm. Gynedil is made entirely of medical grade materials and is gamma sterilized.

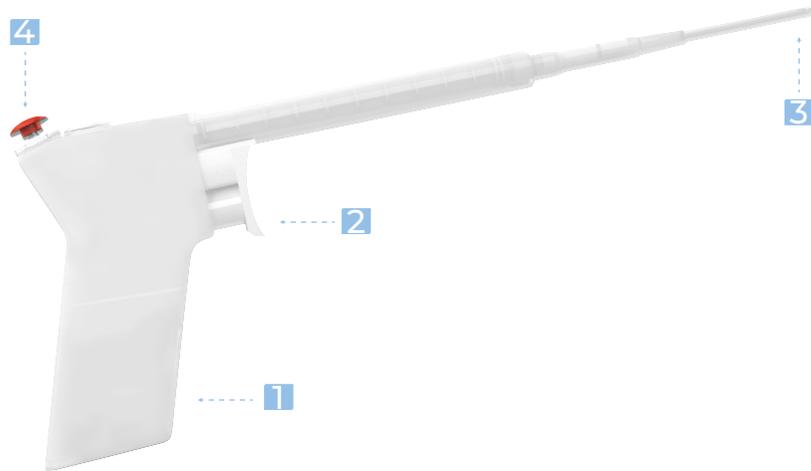

- 1 Reservoir prefilled with distilled water. Gynedil comes ready for immediate use.
- 2 Fast dilation with controlled rhythm and feedback on pressure.
- 3 Flexible dilator tip in medical grade silicon, latex free, causes no allergies or irritations. It can be adjusted to individual shapes of the cervical canal and allows easy insertion.
- 4 Release valve for fast deflation and easy removal of the balloon dilator.

## Gynedil - How it works:

- 1 A flexible dilation tip is inserted into the cervical canal along with its entire length following the cervix anatomy.
- 2 By pressing the trigger (20-22x), the balloon dilator inflates laterally without friction.
- 3 After 30 seconds, dilation is complete. By pressing the upper button on the handle the liquid returns to the handle.
- 4 The balloon dilator deflates to its original diameter and is removed from cervical canal.



Gynedil has a safety system that regulates the pressure, so after obtaining pressure necessary for dilation, it does not allow any further pressure increase, no matter how many times you pull the trigger.